

Specialists in endoscope reprocessing
for more than 30 years!

WASSENBURG® WD440 PT

Pass Through Endoscope
washer-disinfector



- + **NEW LED process identification**
Clear process identification through colour LED technology.
- + **NEW Intuitive colour touchscreen**
User-friendly: fast, safe and easy operation.
- + **NEW High-level ergonomics**
Hygienic one-touch opening by foot sensor.
- + **Outstanding controllability**
Clever clean-contaminated separation with incredibly low air leakage between the clean and contaminated area.
- + **Full traceability**
Access to processes for monitoring, analysis and calculation has never been so easy and reliable by using a data printout or Process Manager.

OPTIMUM CONTROL

- + Channel separation
- + Flow and pressure control
- + Tube connection control
- + Independent Monitoring System (IMS)
- + Leakage control
- + Blockage control
- + Traceability with Process Manager
(and upload possibility to third party programmes)

DIMENSIONS AND WEIGHT


Width x Depth x Height: 1100 x 790 x 2060 mm (+ 15 mm on top, for several connection plugs)
Weight when empty (approx.): 325 kg

ELECTRICAL CONNECTION

Voltage: 220 – 240 V (frequency 50Hz)
Type of electrical connection: plug with ground connection

WATER SUPPLY AND DRAIN

Min. / max. pressure: 2.5 / 4 bar (max. difference between hot and cold water supplies: 1 bar)
Minimum flow rate for hot / cold water supplies: 10 l/min (drain capacity: 20 l/min per drain)
Drain: 2 drains + 2 overflows, all tubes will be shared in 1 drain (diameter drain: >DN50)
Height of drain: 350 mm (+/- 5 mm)

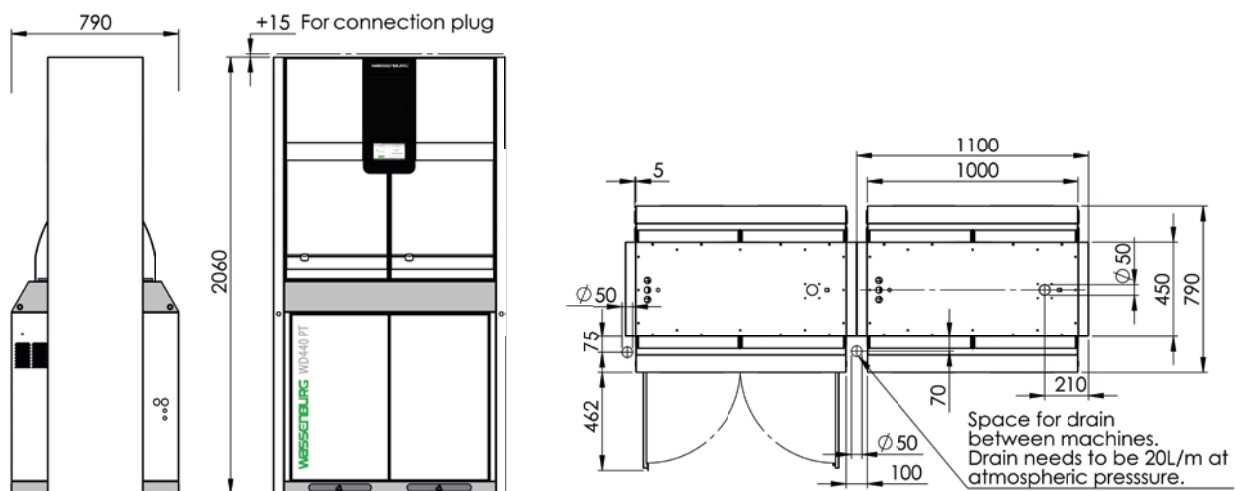

Please contact us for the full installation specifications of the WASSENBURG® WD440 PT.
The WASSENBURG® WD440 PT is CE-marked and fully compliant with EN ISO 15883 parts 1 and 4.
© Copyright 2015. Subject to alterations.