

Specialists in endoscope reprocessing for more than 35 years!




ENDOCART

Transport solution for all types
of flexible endoscopes

A SAFE ENDOSCOPE FOR EVERY PATIENT

Wassenburg Medical is a leader in endoscope reprocessing. For more than 35 years we have been developing innovative solutions with one goal: to provide a safe endoscope for every patient. This has resulted in a range of products and solutions that cover the full spectrum of endoscope reprocessing: a total solution.

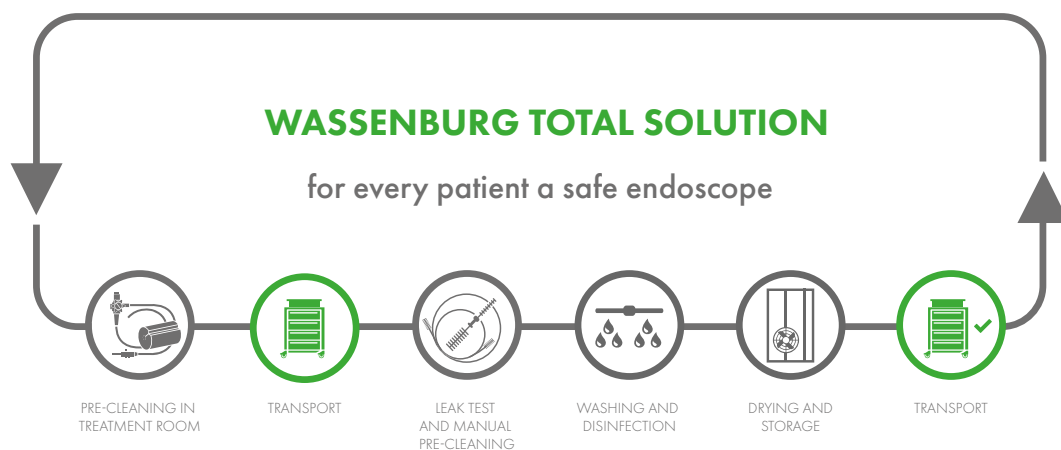

BENEFITS OF THE ENDOCART

- Protecting valuable equipment during transport
- Minimizing the risk of cross contamination
- Protecting medical staff from microbiological hazards
- Improving the workflow in daily operation
- Easy to use; status identification at a glance
- Easy to clean; easy maintenance



- Smooth surface of the tray for ease of cleaning and disinfection
- The top of the trolley provides additional workspace in smaller areas (e.g. patient bedside)
- Tray has no hidden areas for optimum cleaning and disinfection
- No unnecessary coiling of the endoscope because of the tray size
- Trays and top are heat resistant up to 93°C resulting in A0 600 process level according to ISO 15883 part 1



- Removable top for easy cleaning and disinfection
- 6 trays to transport multiple endoscopes from reprocessing to the next patient
- Frame made of stainless steel
- User friendly handle ensures a firm grip pushing the trolley
- A range of two different liners (sterile/non sterile) designed for the EndoCart will protect the disinfected endoscope during transport
- Easy to manoeuvre with four rotating castors
- Equipped with brakes to improve stability



PRODUCT SPECIFICATIONS ENDOCART

The Wassenburg EndoCart is designed to protect flexible endoscopes from airborne contaminants and potential damage during transport between the examination room and the reprocessing facility. The EndoCart robust yet elegant outline in combination with four rotating castors makes it easy to manoeuvre in smaller areas.

Trolley dimensions	510 x 690 x 1050 mm
Capacity	6 trays
Weight (empty)	23 kg
Weight (including trays)	31 kg

PRODUCT SPECIFICATIONS WASSENBURG LINERS AND TRAY

The liner package consists of two latex free liners; a base liner to be placed in the EndoCart Tray to transport reprocessed endoscopes on, and a liner with clear colour coded visual identification to cover the tray to offer a protective environment for both patients as medical staff. The size of the tray prevents endoscopes from kinking and or physical shocks. The tray has no lip portions that could interfere with an effective cleaning. The Wassenburg Tray Liner folds around the tray edges.

Tray dimensions	430 x 550 x 100 mm
Capacity	1 endoscope
Weight (empty)	1.3 kg
Wassenburg Tray Liners	200 items (4 X 50 per box)
Wassenburg Tray Liners - non-sterile	300 items (per box)

Wassenburg Ltd.
1 Smithy Wood Drive
Smithy Wood Business Park
S35 1QN SHEFFIELD
UNITED KINGDOM

T: +44 (0)114 245 4282
sales@wassenburgmedical.co.uk
www.wassenburgmedical.co.uk

© Copyright 2023. Subject to alterations.